

Smoke shark™

Smoke Evacuator

A Sleek, Powerful, Smoke-Eating Machine!



- Easy To Operate
- Quiet and Effective Smoke Evacuation
- Long-life 18-hour Filter
- Low Cost Per Procedure
- Excellent Particle Capture
- Quality Made in the U.S.A.



The Aaron® Smoke Shark™ is designed to be easy to use and extremely effective at smoke and particle removal at the operative sight. While designed for use with all electrosurgical generators, it is also excellent for use in controlling laser smoke. The Smoke Shark's 18-hour filter combines four-stage filtration into one casing (Pre-filter, ULPA, Carbon, Post-filter), keeping cost per procedure down and minimizing waste. The four-stage filter captures particles down to 0.01 μm with 99.99995% efficiency giving you the confidence you need for a safe work environment. The 24 ounces of carbon contained in the Smoke Shark filter; over twenty times the amount of a competing brand, does a superior job of odor elimination.

A variety of accessories are available to customize the Smoke Shark to your specific application. These include laser-resistant tubing, tubing reducers for attachment to smoke evacuating pencils or vaginal speculum, many other tubing configurations (sterile or non-sterile), a foot control, and a remote switch activator.

Smoke Shark Accessories



SF18
long-life 18-hour filter - 1/box



SEAS
arm stand - 1/box



786T
7/8" x 6' tube, non-sterile - 24/box




SETW
7/8" x 6' tube with wand & tip, non-sterile - 24/box



786TS
7/8" x 6' tube, sterile - 24/box



SETWS
7/8" x 6' tube with wand & tip, sterile - 24/box

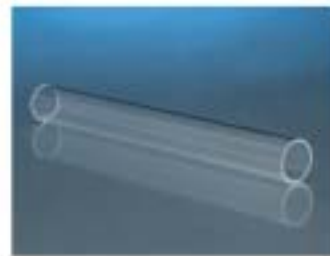

SELW
laser-resistant wand
7/8" x 8", non-sterile - 3/box



SELM
laser mask - 50/box



SELWS
laser-resistant wand
7/8" x 8", sterile - 10/box



SEFM
particle respirator cone mask - 20/box



FSSE
footswitch - 1/box



SERS
remote switch activator - 1/box



SEFC
fluid collection canister 1200cc - 6/box



SELVS
laparoscopic tubing kit with valve, sterile - 12/box



SERH
canister ring holder - 1/box



SERF
reducer fitting 7/8" to 1/4" x 24" - 10/box [fits optional vaginal speculum or Aaron® ESP4 smoke & fluid evacuating pencil]



SEWP
wall plate to hold canister - 1/box



SERFS
reducer fitting, sterile - 10/box



SELKS
laparoscopic tubing kit, sterile - 12/box



SEVL
7/8" x 6' tubing reduced to 3/8" for ValleyLab & IEC pencil, sterile - 5/box

STANDARD OFFICE SETUPS

ELECTROSURGERY

1 - SE01 Smoke Shark™ Unit
 1 - SF18 Long-life 18-hour Filter
 1 Case (24) - 786T 7/8" x 6' Tube*
 or 1 Case (24) - SETW 7/8" x 6'
 Tube w/Wand*
 1 - FSSE Footswitch
 1 - SEAS Arm Stand

OBGYN

1 - SE01 Smoke Shark™ Unit
 1 - SF18 Long-life 18-hour Filter
 1 Case (24) - 786T 7/8" x 6' Tube*
 or 1 Case (24) - SETW 7/8" x 6'
 Tube w/Wand*
 1 - FSSE Footswitch
 1 Case (10) - SERF Reducer Fitting
 7/8" to 1/4" x 24"

LASER

1 - SE01 Smoke Shark™ Unit
 1 - SF18 Long-life 18-hour Filter
 1 Case (24) - 786T 7/8" x 6' Tube*
 1 Box (3) - SELW 7/8" x 8"
 Laser Resistant Wand*
 1 - FSSE Footswitch
 1 - SEAS Arm Stand
 *Non-Sterile

Two Complete Systems

The Aaron Smoke Shark is available with two complete systems specifically designed for the OBGYN:



The Aaron 1200-G

The Aaron 1200™ is a multipurpose electrocautery generator for use in the modern physician's office and surgical center. It features five operating modes and 120 watts of cutting power. The autoclavable single-button handpiece holds all 2.3 mm (3/32") shaft electrodes securely for the greatest flexibility and dependability available.

A1200-G INCLUDES:

- **A1200 Electrocautery Generator**
 - Five Extra Disposable Electrodes
 - A1202C Reusable Grounding Cord
 - 1202 Five Solid Grounding Pads
 - A801 Reusable Single Button Pencil
- **SEGYN Smoke Shark System**
 - SE01 Smoke Shark
 - SF18 Eighteen Hour Filter
 - 786T 7/8" Tube 24" Long
 - FSSE Footswitch
 - SERF Reducer Fitting (Attaches to Vaginal Speculum)
- **ESLK Electrocautery Loop Kit**
 - ES09 10mm x 10mm Loop
 - ES11 15mm x 10mm Loop
 - ES12 20mm x 8mm Loop
 - ES13 20mm x 15mm Loop
 - ES16 10mm x 8mm Square
- **A812-C Electrocautery Mobile Stand with Trays and Height Extender**



The Aaron 1250-G

The Aaron 1250™ is a multipurpose electrocautery generator incorporating automatic safety features like self-test circuits, audible tones, discreet outputs, isolated circuitry, and Bovie NEM™ (neutral electrode monitoring) pad sensing. It features 120 watts of cutting power, five modes of operation, and digital error detection to satisfy all surgical demands with safety, reliability, and convenience.

A1250-G INCLUDES:

- **A1250 Electrocautery Generator**
 - Five Extra Disposable Electrodes
 - Two Disposable Pencils
 - A1252C Reusable Grounding Cord
 - ESRE-1 Five Split Grounding Pads
- **SEGYN Smoke Shark System**
 - SE01 Smoke Shark
 - SF18 Eighteen Hour Filter
 - 786T 7/8" Tube 24" Long
 - FSSE Footswitch
 - SERF Reducer Fitting (Attaches to Vaginal Speculum)
- **ESLK Electrocautery Loop Kit**
 - ES09 10mm x 10mm Loop
 - ES11 15mm x 10mm Loop
 - ES12 20mm x 8mm Loop
 - ES13 20mm x 15mm Loop
 - ES16 10mm x 8mm Square
- **ESMS-C Electrocautery Mobile Stand with Smoke Evacuator Tray and Instrument Drawer**

Distributed By:



7100 30th Avenue North
St. Petersburg, FL 33710-2902
U.S. Phone 1-800-537-2790 • Fax 1-800-323-1640
Int'l Phone +1-727-384-2323 • Fax +1-727-347-9144
www.aaronmed.com • sales@aaronmed.com